


Customer Information Aid

Double-Donning

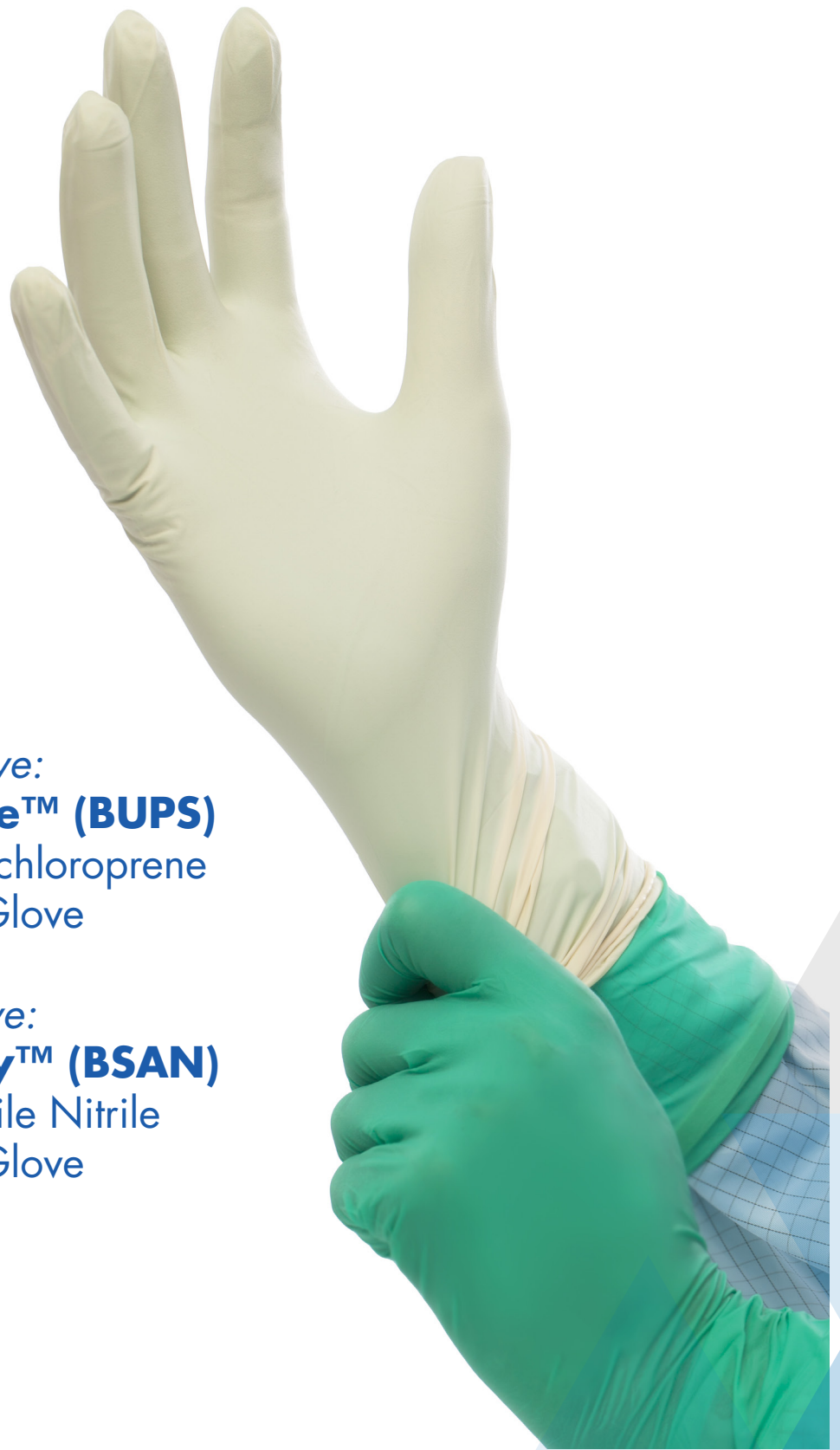


Outer Glove:
BioClean Ultimate™ (BUPS)
300mm Sterile Polychloroprene
Cleanroom Glove

Inner Glove:
BioClean Synergy™ (BSAN)
300mm Non-sterile Nitrile
Cleanroom Glove

Europe (Headquarters)
Americas
Asia & Middle East

europa@bioclean.com
americas@bioclean.com
malaysia@bioclean.com

www.bioclean.com

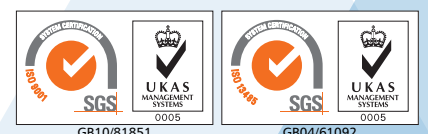

Customer Information Aid

Double-Donning



- 1** Don a pair of BioClean Synergy™ (BSAN) non-sterile powder-free gloves. Ensure the cuffs are pulled up to their maximum length.



- 2** Take a pair of BioClean Ultimate™ (BUPS) sterile powder-free gloves and tear open the 'EasyTear™' pouch at the top where indicated.



- 3** Remove the inner plastic wallet, place on a flat surface and unfold so that both gloves are accessible. Do not touch the exterior of the gloves. Pinch the sides of the plastic wallet to ensure flaps do not hamper access to the gloves.

- 4** Insert fingers into the first BUPS glove, palm side up, bend thumb towards the palm and slide hand fully into glove whilst pulling up on the cuff, ensuring no contact with the outside surface of the glove (i.e. non-sterile against non-sterile contact). Leave the cuff folded at this stage.



- 5** Place fingers of donned glove hand into cuff of BUPS glove left in inner plastic wallet to give stability. Slide fingers of BSAN gloved hand into BUPS glove, bend thumb towards palm and slide fingers fully into glove whilst pulling glove up over the wrist using gloved hand touching the outside surface only (i.e. sterile against sterile contact).

- 6** Pull up the cuff of the first BUPS glove using the BUPS gloved second hand. Take care to only touch the outside of the glove. (i.e. sterile against sterile contact).